



PERFECTCOLL™

OPERATING MANUAL

# PC FACE GOLD 24 K

Soft Skin Effect





# Contents



24 carat gold of the highest hallmark	4
Recommended cosmetics	6
Design and technical data	10
Application	12
Rejuvenation effects	14
Indications and contraindications	18
Care and maintenance	19

# 24 carat gold of the highest hallmark




**PC FACE GOLD 24 K**, made with the latest technology, is equipped with a head covered with 24 carat gold. **PC FACE GOLD 24 K** has a specially designed shape and vibrates **6 000 times** per minute. The device releases the face, neck and décolleté muscle tension. In addition, it perfectly improves firmness and elasticity of flabby skin.

It is equipped with a vibration system, which stimulates the deeper layers of the skin, helping to improve its blood supply, firmness and flexibility.



In this way, there is a phenomenon of increased proliferation of fibroblasts. Thus, new collagen and elastin fibres are produced, which provide the skin with firmness and elasticity. As a result, the skin is younger and healthier. These processes are the result of the “regenerative properties” of gold. Pure gold ions are also referred to as active gold. Gold ions can be absorbed through the skin with the cosmetic.

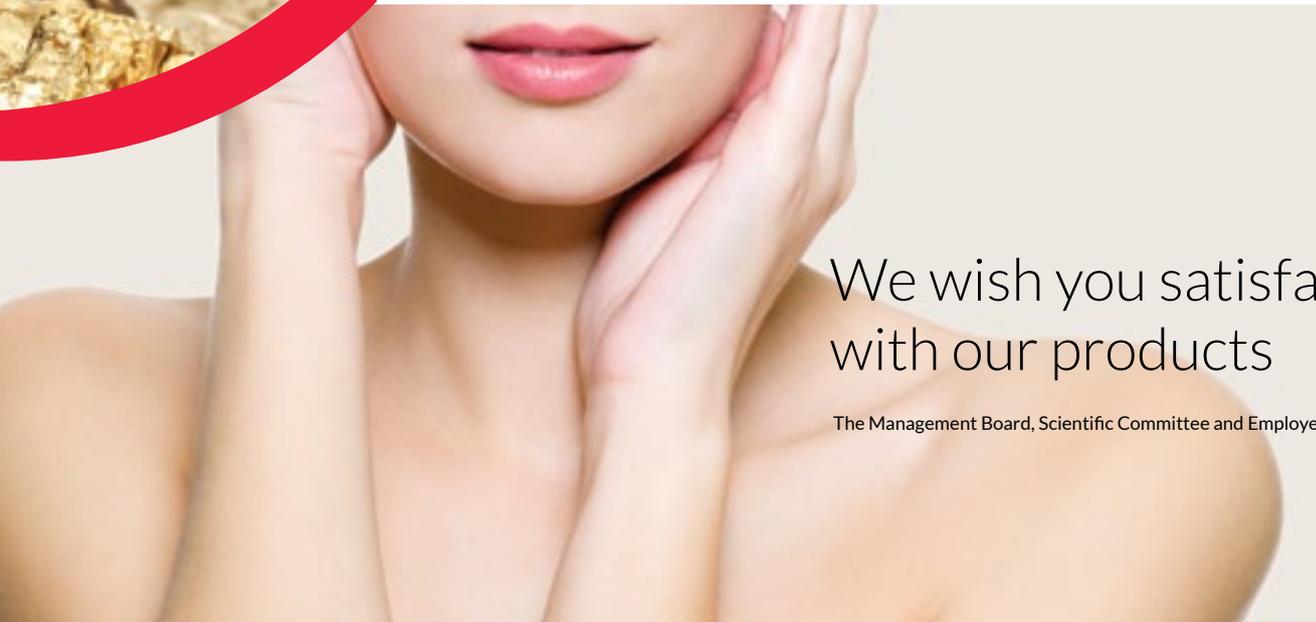

We wish you satisfaction  
with our products

The Management Board, Scientific Committee and Employees of Perfect Coll SP. z o.o.

# Recommended cosmetics

Use **PC FACE GOLD 24 K** only with the original cosmetics of PERFECT COLL



1

Before treatment

## PC DETOX EXTRACT FROM NORWEGIAN WILD CORNFLOWER

The cleanest and most modern, all-natural extract of the flowers of Centaurea harvested from organic farming in Scandinavia, enriched with natural enzymes from Organic Kelp (organic algae). It gently cleans and soothes the sensitive skin of face, neck and the whole body. It does not contain parabens, chemical preservatives, and drying skin alcohol.

## SYMBIOSIS OF SCIENCE WITH NATURE

PERFECT COLL cosmetics have unique, innovative formula developed in biotechnology, environmental laboratories in the Polish medical universities.

We are keen on nature and ecology, so the cosmetics with the PC sign are healthy for your skin.



2

During treatment



### PC PREBIOTIC HYALURONIC SERUM

It has innovative and entirely natural prebiotic formula and the most efficient ratio of the various particle sizes of hyaluronic acid, which causes deep hydration of dehydrated or desiccated skin layers. Due to the unique content of prebiotic enzymes, serum is a perfect way for the most demanding skin. It is the only prebiotic product on the market characterized with extremely high absorption.



3 After treatment

### PC PLATINUM COLLAGEN +32°C HIALURON


Poland is famous for producing the purest and best in the world collagen perfectly assimilated by the human body. Our collagen, thanks to the unique technology, has the unique heat-resistance to +32°C. As the only product on the market it may hold fully qualified name of the Platinum, because it contains platinum biopeptide. Platinum facilitates penetration of biologically active substances into the deep layers of the skin, optimizing action of the hyaluronic acid present in the preparation. Excellent composition of three active components in hyaluronic acid enhances the perfect effect after treatment with **PC FACE GOLD 24 K**.



4 After treatment

### PC PREBIOTIC SOOTHING CREAM

Natural and eco-friendly, innovative soothing cream, which due to unique composition has remarkable properties to support the work of damaged skin cells. Its components are responsible for regeneration of even the most delicate skin, and enzymes from papaya fruit provide the effect of light glowing. It is the only product on the market that contains prebiotic enzymes, which provide very high absorbency.

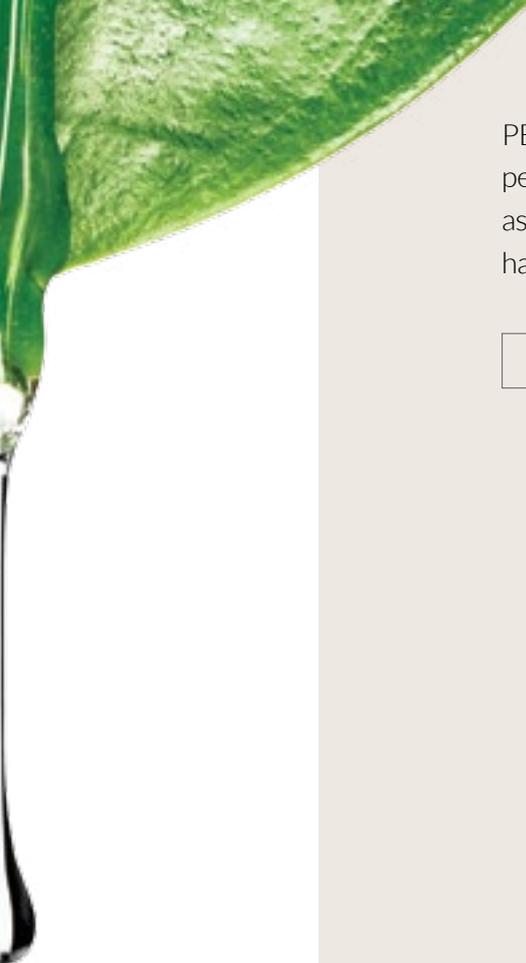

PERFECT COLL cosmetics do not contain parabens, chemical preservatives, petroleum (paraffin), mineral oils, alcohols, and glycols drying out the skin, as well as irritating foaming agent Sodium Laureth Sulfate and other chemical substances harmful to the skin.

Refer to "Perfect Coll Products Catalogue"



# Design and Technical data PC FACE GOLD 24 K





## TECHNICAL DATA

### DEVICE:

Name of the product	PC FACE GOLD 24 K
Dimensions	140 x 40 mm
Weight	81 g
Power supply	1 x AA

MADE IN PRC FOR PERFECT COLL



# Application

# PC FACE GOLD

# 24 K



## SIMPLE STEPS TO USE PC FACE GOLD 24 K:

- Use the AA battery and place it properly in the handle of the device. To turn on the device, use the tip and turn it to ON. To turn off the device, turn the tip to OFF.
- Remove make-up.
- Wash the face with **PC DETOX EXTRACT FROM NORWEGIAN WILD CORNFLOWER**.
- Apply the right amount of **PC PREBIOTIC HYALURONIC SERUM**.
- Then hold the gold head on the face, gently and slowly move from the middle of the chin to outside. If there are wrinkles, you can massage longer, for example, in the vicinity of eyes, mouth, or neck.
- After the treatment, wash your face with **PC DETOX EXTRACT FROM NORWEGIAN WILD CORNFLOWER**.
- Apply **PC PLATINUM COLLAGEN +32°C HIALURON**, and massage.
- After the procedure, use **PC PREBIOTIC SOOTHING CREAM**.
- Clean the gold spoon (*Chapter Care and maintenance*).
- Put **PC FACE GOLD 24 K** into the box.



### NOTE!

Keep the device at 1/3 of the height of the handle. If you hold too high, this can result in release of vibration, and thus weaken the effect of the treatment. In hard-to-reach places, for example, around the nose wings, you can use the rounded edges of the head. You can apply spot massage to the chosen areas of the body. Hold the gold head in one place for approx. 2-3 seconds.



Frequency of treatment:  
Every day




Treatment time:  
10 min.



# Effects of rejuvenation

- Lymphatic detox
- Lifting
- Reduced wrinkles
- Slim face and improvement of its contours
- Lightening and smoothing
- Rejuvenating and anti-ageing effects

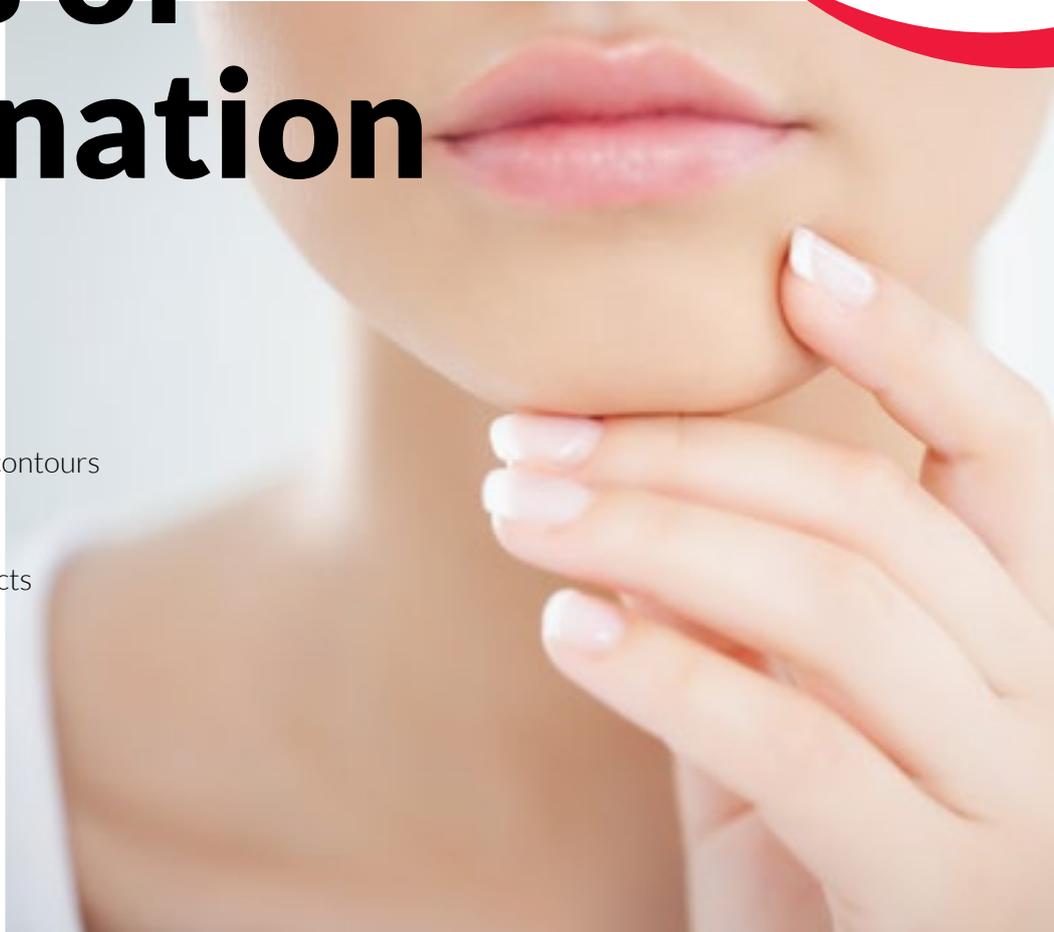



In the human body there is poor electrical conductivity.

In healthy and young skin, current normally flows through the vessels in the blood and lymphatic circulations, improving metabolism, so your face is supple and full of vitality. Unfortunately, as we age, changes in the body cause ion imbalance, which has negative consequences in the form of bold stratum corneum of the epidermis, microcirculation loss and blocked electrical conductivity. Metabolic wastes and water excess cannot be discharged and are accumulated in the skin. Consequently, cells activity decreases and the skin problems appear, such as swelling of the face, wrinkles, dehydration, skin sagging and discolourations.

A unique T-shaped head of **PC FACE GOLD 24 K** is covered with 24 carat gold. In addition, it is equipped with the system of vibrations, frequency of which is **6000 times per minute**. The vibrations stimulate the skin, respectively, because the negative ions of gold caused by vibrations have the same wavelength, as electricity in the human body. For this reason, current flow in your body will improve, thanks to a gentle stimulation by **PC FACE GOLD 24 K**. In summary, the potential of the electrical balance stimulates cell metabolism, reduces swelling of the skin, activates the cells to renew, lightens pigmentation, reduces wrinkles. It makes the skin not only lighter, but also moisturised, healthy, tense and visibly rejuvenates it.



## LYMPHATIC DETOX


The head of **PC FACE GOLD 24 K** is covered with 24 carat gold. Its unique shape allows you to execute the vibrations stimulating microcirculation and skin cells. It balances energy potential and improves lymphatic circulation in the skin. In addition, it renews functions of deeper cells, structural fibres, causes oxygenation of the cells, exchange of ions and increases absorption of nutrients. It efficiently removes unnecessary products of cell transformations and toxins. As a result, your skin regains vitality and youthful glow.

## LIFTING EFFECT

Electrical voltage in the normal cells is in the range from 0.06 to 0.07 V. Taking into account the negative impact of the environment, air pollution, or accompanying daily stress, the voltage drops from 0.02 to 0.03 V. As a result, the skin loses its tightness and elasticity. **PC FACE GOLD 24 K**, due to its unique vibration feature, stimulates cells activity, improves their energy, which causes lifting.

## REDUCING WRINKLES

**PC FACE GOLD 24 K** emits negative ions, which stimulates cells in the deeper layers of the skin. Thanks to the unique features, **PC FACE GOLD 24 K** strengthens collagen fibres, the skin layers become thicker, which increases skin flexibility. All these processes in an unusual way contribute to the shallowing of wrinkles and to improve the appearance of the face.

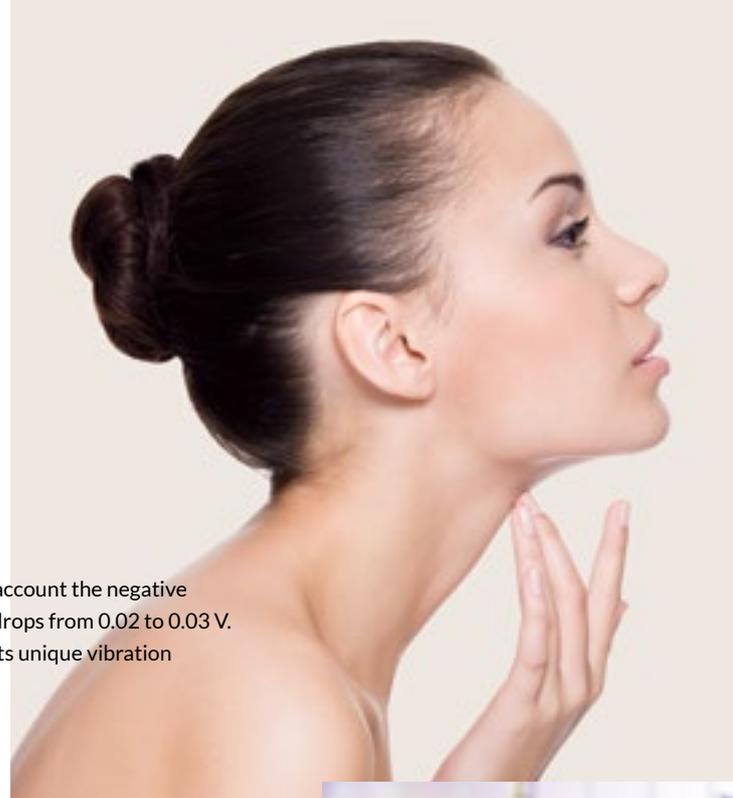

## SLIMMING THE FACE AND IMPROVING ITS CONTOURS

**PC FACE GOLD 24 K** emits **6000 vibrations** per minute. Massage stimulates blood microcirculation, increasing the temperature of the deeper layers of the skin. This definitely increases the activity of cells, improves muscle tone and removes the excess of body fat.

## LIGHTENING AND SMOOTHING

T-shaped **PC FACE GOLD 24 K** improves microcirculation and activates skin cells. It increases protein synthesis and oxygenation of the skin, enhances exchange of ions and enables deeper penetration of nutrients. In addition, it helps to eliminate unnecessary cell transformation products and toxins. It improves microcirculation on your face and makes your skin glow and smoother.

## REJUVENATING AND ANTI-AGEING EFFECT

The head of **PC FACE GOLD 24 K** emits negative ions flow. Negative ions neutralize the effects of free positive ions, thus blocking the effects of free radicals. It increases blood circulation, improves metabolism, making your skin smooth, stimulates it and provides the effect of visible rejuvenation.



# Indications

- Lymphatic stasis.
- Ageing skin.
- Skin with tendency to contamination.
- Swells and dark circles under the eyes.
- Wrinkles around the eyes, and, in particular, "crow's feet".
- Skin with discolouration.
- Lack of firmness and elasticity of the skin.
- Dry, grey skin, without glow.



# Contraindications

- Sun-burned skin (in the place of treatment).
  - Break in the epidermis (in the place of treatment).
  - Bacterial, viral, fungal infection (in the place of treatment).
  - Fresh scars, birthmarks and tattoos (in the place of treatments).
  - Severe eczema, psoriasis, atopic skin, rosacea (active lesions in the place of treatment).
  - Allergy to the components of the cosmetics.
  - Surgery provided within the last 6 months (in the place of treatment).
  - Treatments including administration of botulinum toxin provided within the last 9 months (e.g. Botox) (in the place of treatment).
  - Corticosteroids (in the place of treatment). Autoimmune diseases (e.g. HIV, lupus, etc.).
  - Cancer.
  - Tumours, nodules of unknown origin.
- Pregnancy.
  - Lactation (breastfeeding).
  - Increased temperature.
  - Cardiovascular disease.
  - Epilepsy.
  - Tumours, nodules of unknown origin.
  - Arterial hypertension.
  - Tics.
  - Cardiac dysrhythmia.

# Care and maintenance

- After the treatment, remove the remaining preparation from the head.
- Wash **PC FACE GOLD 24 K** under running water.
- Wipe with a tissue soaked with antibacterial preparation (e.g. 2% salicylic spirit).
- Put the device into the box.
- If you do not intend to use the device for an extended period of time, take out the battery.



Keep the device clean and dry.





## Centrum Operacyjno-Rozwojowe Perfect Coll Sp. z o.o.

Pomorski Park Naukowo-Technologiczny (lok. A 210)  
Al. Zwycięstwa 96/98, 81-451 Gdynia



[kontakt@perfect-coll.pl](mailto:kontakt@perfect-coll.pl)



+48 58 767 42 24



[www.perfect-coll.pl](http://www.perfect-coll.pl)

